



Shanghai Genomics, Inc.

No. 647 Songtao Road

Shanghai 201203, China

Phone: 86-21-50802786

FAX: 86-21-50802783

www.shanghaigenomics.com

TOTAL RNA- HUMAN BLADDER

CATALOG NUMBER SG1110-11

QUANTITY 50 µg

DESCRIPTION Total RNA from normal human bladder tissue.

DONOR INFORMATION Gender: Male Age: 35 Race: Asian Cause of Death: Acute Head Injury

QUALITY CONTROL The OD 260/280 ratio for total RNA sample in Tris-HCl buffer (pH 7.5) is 2.0. RNA integrity is confirmed by the electrophoresis of 1 µg total RNA sample on a denaturing formaldehyde/agarose/EtBr gel (Fig.1).

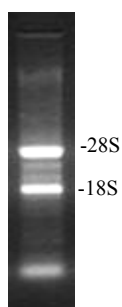


Figure 1: Total RNA from human bladder. 1 µg of total RNA was analyzed using a 1% denaturing formaldehyde /agarose /EtBr gel.

STORAGE BUFFER 50µg RNA/vial is dissolved in 70% ethanol, 0.08M sodium acetate.

PRIOR TO USE

- 1) Centrifuge the tube at 14,000 rpm for 10 minutes in a microcentrifuge at 4°C
- 2) Remove the supernatant and resuspend the pellet in 70% ethanol.
- 3) Centrifuge the tube at 14,000 rpm for 10 minutes in a microcentrifuge at 4°C.
- 4) Remove the supernatant and resuspend the pellet in DEPC-treated H₂O to the desired concentration.

STORAGE/HANDLING Purified RNA is recommended being stored at -70°C. Avoid repeated freeze /thaw cycles. Handle with gloves to prevent contamination with RNase.

SHELF LIFE 1 year from date of receipt under proper storage conditions.
