



Shanghai Genomics, Inc.

No. 647 Songtao Road

Shanghai 201203, China

Phone: 86-21-50802786

FAX: 86-21-50802783

www.shanghaigenomics.com

## TOTAL RNA- PAIR OF HUMAN ESOPHAGUS TUMOR

---


**CATALOG NUMBER:** SG1130-03

**QUANTITY:** 35 µg x 2

**DESCRIPTION:** Total RNA from pair of human esophagus tumor and normal adjacent tissue.

**DONOR INFORMATION:** Gender: Male Age: 57 Race: Asian Pathology Diagnosis: Esophageal squamous cell carcinoma

**QUALITY CONTROL:** The OD 260/280 ratio for total RNA sample in Tris-HCl buffer (pH 7.5) is 2.0. RNA integrity is confirmed by the electrophoresis of 1 µg total RNA sample on a denaturing formaldehyde/agarose/EtBr gel (Fig.1).

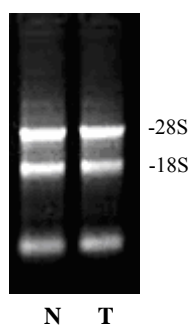

**Figure 1: Total RNA from pair of human esophagus tumor and normal adjacent tissue.** 1 µg per lane of total RNA was analyzed using a 1% denaturing formaldehyde /agarose /EtBr gel.

**STORAGE BUFFER:** 35 µg RNA/vial is dissolved in 70% ethanol, 0.08M sodium acetate.

**PRIOR TO USE:**

- 1) Centrifuge the tube at 14,000 rpm for 10 minutes in a microcentrifuge at 4 °C
- 2) Remove the supernatant and resuspend the pellet in 70% ethanol.
- 3) Centrifuge the tube at 14,000 rpm for 10 minutes in a microcentrifuge at 4 °C.
- 4) Remove the supernatant and resuspend the pellet in DEPC-treated H<sub>2</sub>O to the desired concentration.

**STORAGE/HANDLING:** Purified RNA is recommended being stored at -70 °C. Avoid repeated freeze /thaw cycles. Handle with gloves to prevent contamination with RNase.

**SHELF LIFE:** 1 year from date of receipt under proper storage conditions.

---