



Shanghai Genomics, Inc.  
No. 647 Songtao Road  
Shanghai 201203, China  
Phone: 86-21-50802786  
FAX: 86-21-50802783  
www.shanghaigenomics.com

## *VprBP Polyclonal Antibody*

---

**CATALOG NUMBER:** SG4220-28  
**QUANTITY:** 100µl  
**SOURCE:** Rabbit  
**DESCRIPTION:** Polyclonal antibodies are produced by immunizing rabbits with a peptide corresponding to amino acid residue 1493-1507 (DNSDLEDDIILSLNE) of human VprBP coupled to KLH.  
**SPECIFICITY/SENSITIVITY:** VprBP antibody detects overexpressed human VprBP (Fig.1).

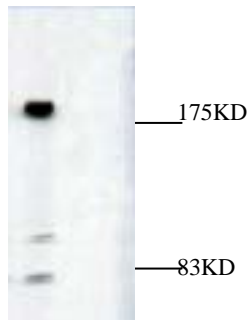


Figure 1: Western Blot analysis of cell lysate expressing VprBP.

**APPLICATION:** Western Blot  
**FORMAT:** Rabbit serum containing 0.02% sodium azide.  
**RECOMMENDED ANTIBODY DILUTION:** Western blot: 1:5000.  
**STORAGE/HANDLING:** Antibody is recommended being stored at -20 . Avoid repeated freeze-thaw cycles.  
**USAGE:** FOR RESEARCH USE ONLY. NOT FOR HUMAN USE.